




New device generation classic-cart

NEW PRODUCT

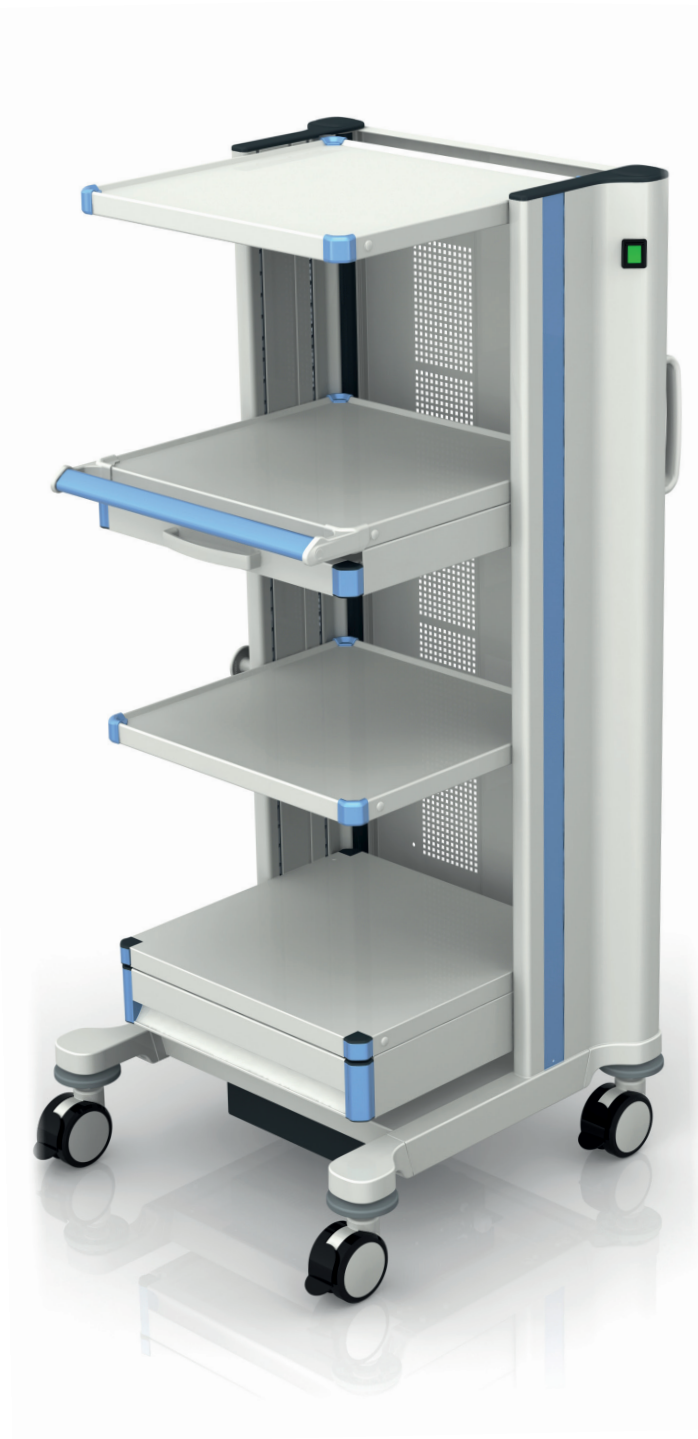

TECHNICAL DETAILS

Load capacity

Total load capacity	150 kg / 330 lbs
Load capacity of shelf	50 kg / 110 lbs distributed load

Dimensions - Basic frame

Width in mm	663 / 763 / 1263
Depth in mm	660
Total height in mm	840 / 968 / 1128 / 1288 / 1448 / 1608 / 1768

Dimensions - Shelf

Width in mm	500 / 600 / 1100
Depth in mm	491

Electrification

Isolating transformer	1200 VA / 1600 VA / 2000 VA
Central on / off switch	yes
Socket strips	yes

Get in touch >>

Europe: +49 89 614425-0 | sales@itd-cart.com
 America: +1 800 947-3901 | salesusa@itd-cart.com
 China: +86 21 5031 8818 | saleschina@itd-cart.com
 India: +91 86 0015 4064 | salesindia@itd-cart.com
 Australia: +61 3 9330 5050 | salesaustralia@itd-cart.com

www.itd-cart.com




DS.0270.799 - 10/2024. Printed in Germany. All rights to make changes reserved.

CUSTOMIZED DESIGN

New colours

BASE COLOR



RAL 9016
Traffic white

DECORATIVE TRIM COLORS



RAL Design 2507025
Meadow blue



RAL 5017
Traffic blue




RAL 7016
Anthracite grey

New design



PERSONALIZED CORPORATE DESIGN



HIGH QUALITY DESIGN WITH OPTIMIZED FUNCTIONALITY

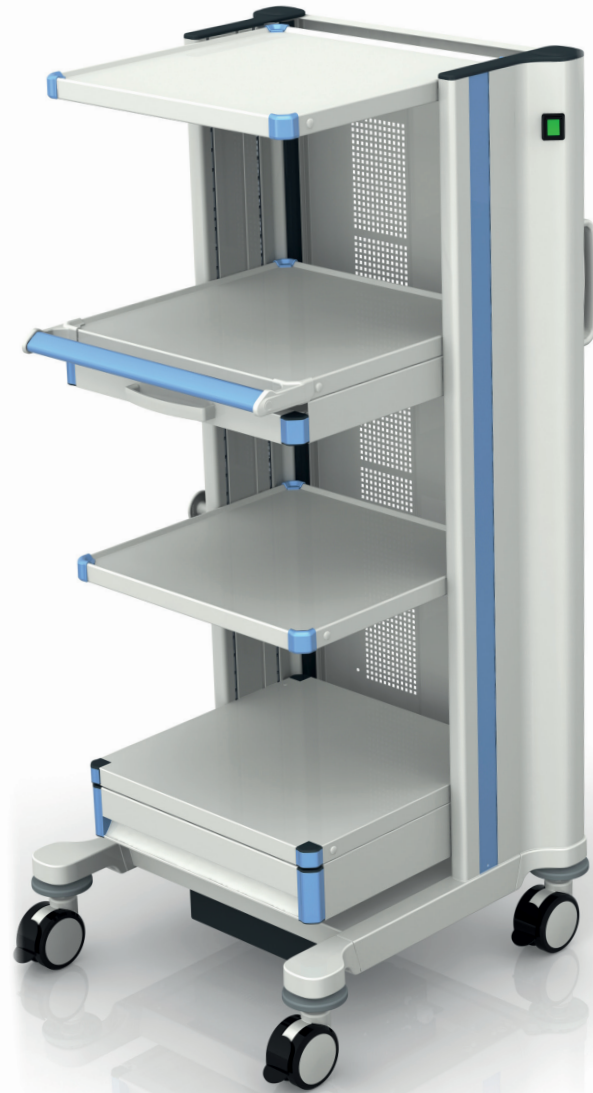


New corners – elegant design with improved hygiene properties

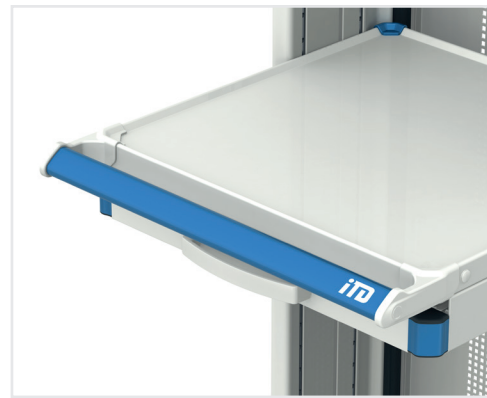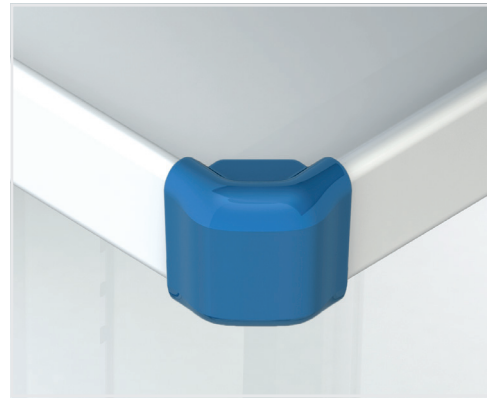

Push handle – comfortable feel, customizable to corporate brand standards

Single, smooth column – offers excellent hygiene features and a distinct design

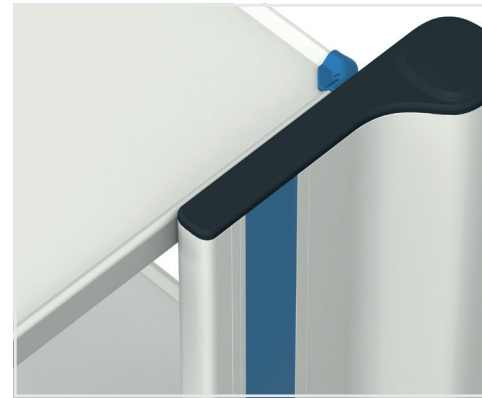

Customized design – highly recognizable





New arm position – large range of motion

Additional rear handles – easy alignment of the cart from any position

Air-permeable rear door – prevents heat congestion of the equipment and unwanted access

International socket strips – unrestricted accessibility

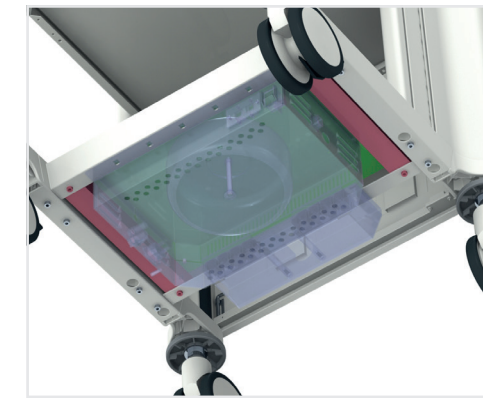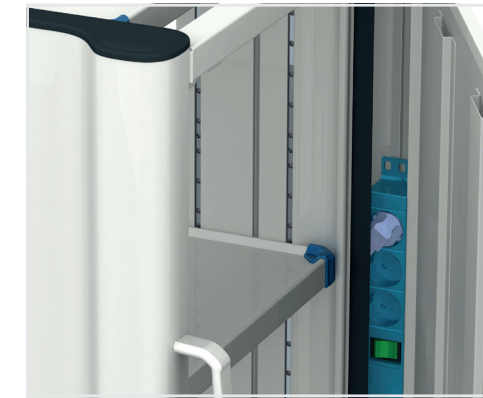

Counterweights – invisibly integrated and flexibly combined with isolating transformers

