



# Medical cart pro-cart

NEW COLORS



## TECHNICAL DETAILS

### Load capacity

Total load capacity	80 kg / 176 lbs
Load capacity of shelf	20 kg / 44 lbs Flächenlast

### Dimensions - Basic frame

Width in mm	610
Depth in mm	620
Total height in mm	884 / 1172 / 1332

### Dimensions - Shelf

Width in mm	450
Depth in mm	400

### Electrification

Isolating transformer	660 VA / 1200 VA / 1600 VA
Central on / off switch	yes
Socket strips	yes

Get in touch »

Europe: +49 89 614425-0 | sales@itd-cart.com  
 America: +1 800 947-3901 | salesusa@itd-cart.com  
 China: +86 21 5031 8818 | saleschina@itd-cart.com  
 India: +91 86 0015 4064 | salesindia@itd-cart.com  
 Australia: +61 3 9330 5050 | salesaustralia@itd-cart.com

www.itd-cart.com



DS-0286.799 - 10/2024. Printed in Germany. All rights to make changes reserved.

## CUSTOMIZED DESIGN

### New colors

#### BASE COLOR



RAL 9016  
Traffic white

#### DECORATIVE TRIM COLORS



RAL Design 2507025  
Meadow blue



RAL 5017  
Traffic blue




RAL 7016  
Anthracite grey

### Color design

For an attractive color design, decorative elements such as handles can be powder-coated in selected colors. The implementation of individual corporate design colors is also possible.



# MOUNTING LARGE MONITORS IN ERGONOMIC POSITION

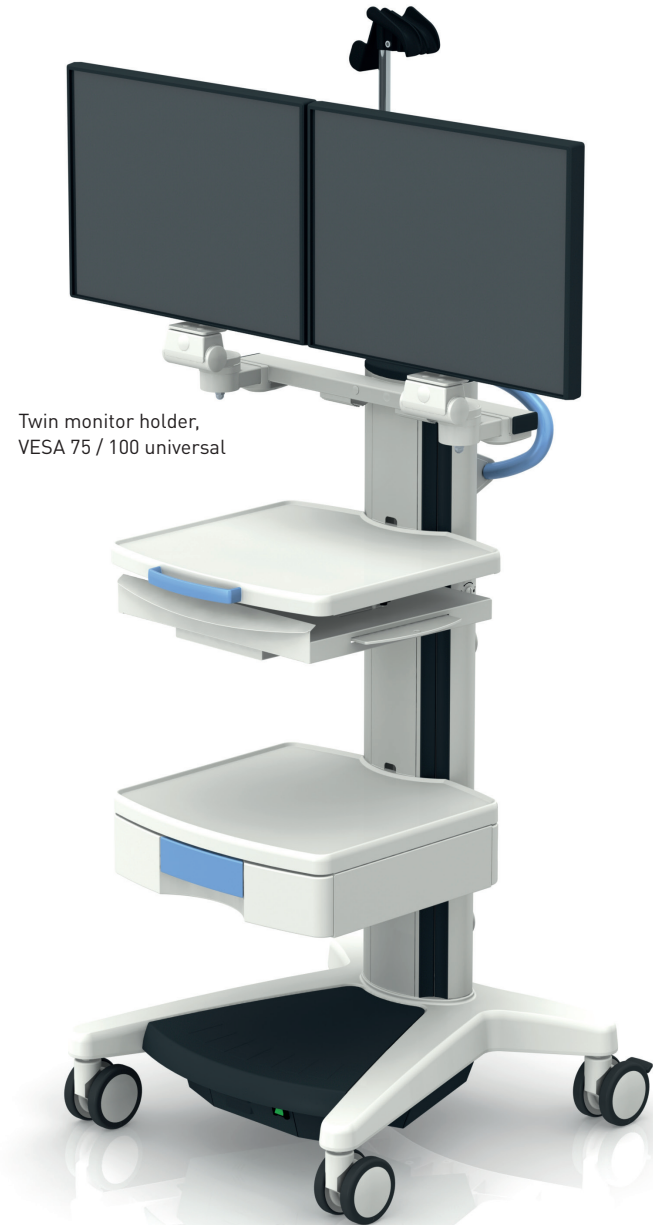

Twin monitor holder,  
VESA 75 / 100 universal

**Flexible monitor holders –**  
for monitors up to 18 kg / 39.6 lbs



Universal adapter

VESA 75 / 100

VESA 75 / 100  
adjustable




**Variable height support arms –**  
high range of movement and  
ergonomic positioning



**The central on/off switch –**  
ergonomic operation of the  
isolating transformer



**Isolating transformer –**  
electrical safety for  
patients and users,  
international use  
through switching the  
voltage

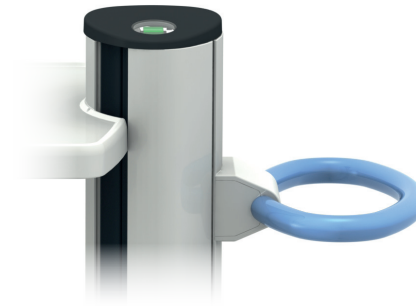

**The foldable push handle –**  
effortless maneuvering and  
space-saving parking position



# SAFE ELECTRIFICATION AND CONCEALED CABLE ROUTING

**Spacious, lockable support column –**  
space for long socket strips, protection  
against unauthorized access and easy  
access to cables



**Flexible rubber strips –**  
cable entry and exit at any position  
along the support column

**Safe cable routing –**  
separate data and power cables  
guarantee electrical safety and  
prevent electromagnetic interference

