



Medical cart vexio-cart

NEW COLORS



TECHNICAL DETAILS

Load capacity

| | |
|------------------------|---------------------------------|
| Total load capacity | 65 kg / 143 lbs |
| Load capacity of shelf | 20 kg / 44 lbs distributed load |

Dimensions - Basic frame

| | |
|--------------------|-------------------------|
| Width in mm | 530 |
| Depth in mm | 610 |
| Total height in mm | 830 / 998 / 1160 / 1320 |

Dimensions - Shelf

| | |
|-------------|-----|
| Width in mm | 424 |
| Depth in mm | 430 |

Electrification

| | |
|-------------------------|--------|
| Isolating transformer | 660 VA |
| Central on / off switch | yes |
| Socket strips | yes |

Get in touch »

Europe: +49 89 614425-0 | sales@itd-cart.com
 America: +1 800 947-3901 | salesusa@itd-cart.com
 China: +86 21 5031 8818 | saleschina@itd-cart.com
 India: +91 86 0015 4064 | salesindia@itd-cart.com
 Australia: +61 3 9330 5050 | salesaustralia@itd-cart.com

www.itd-cart.com



DS-0284.799 - 10/2024. Printed in Germany. All rights to make changes reserved.

CUSTOMIZED DESIGN

New colors

BASE COLOR



RAL 9016
Traffic white

DECORATIVE TRIM COLORS



RAL Design 2507025
Meadow blue



RAL 5017
Traffic blue



RAL 7016
Anthracite grey

Color design


For an attractive color design, decorative elements such as handles and side cable ducts can be powder-coated in selected colors. The implementation of individual corporate design colors is also possible.



ERGONOMICALLY WORKSTATION




The push handle –
enables effortless maneuvering and
can be customized to corporate design
colors

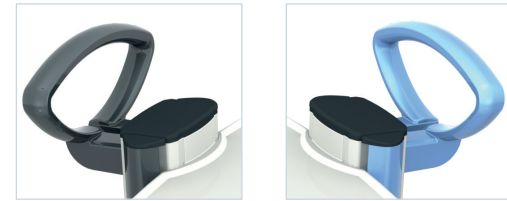

Mechanical height adjustment –
enable a quick change between
standing and sitting workstations

FLEXIBLE ACCESSORY POSITIONING & ORGANIZED CABLE ROUTING

Different column geometries –
enable optimum routing of all cables
along or within the column

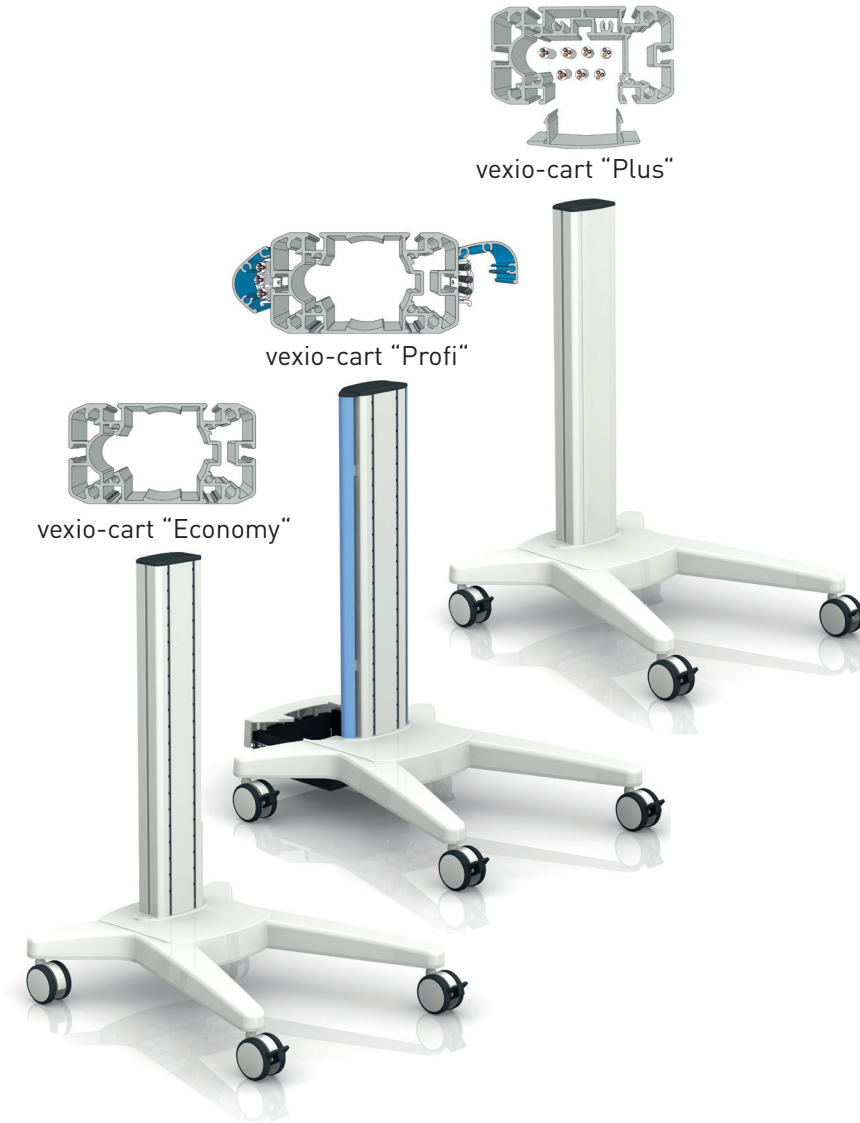


Flexibility when mounting accessories –
the T-slots on all four sides of the profile
offer flexibility in the positioning of
accessories and system components



Isolating transformer 660 VA –
convenient switching on and off at
the end of the column



Global distribution –
various plug types such as cold
appliance, Schuko, French, UK and US
plugs enable global use

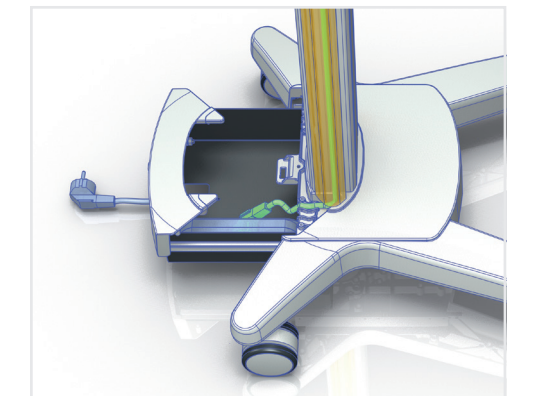

Lockable electrical box –
ensures secure storage of multiple so-
ckets and power supply units and
protects against unauthorized access